



CAPTIVE® J-STRAIGHTENER INSTRUCTIONS FOR USE

TO RELEASE GUIDEWIRE FROM DISPENSER:


1. Grasp the shoulder of the J-Straightener between the forefinger and thumb.

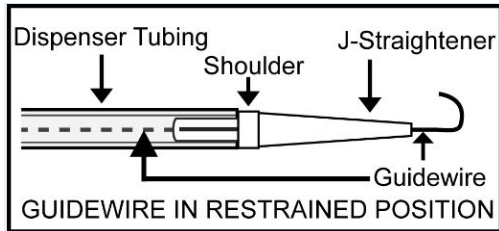

2. Gently pull on the J-Straightener shoulder until it is just removed from the dispenser tubing.

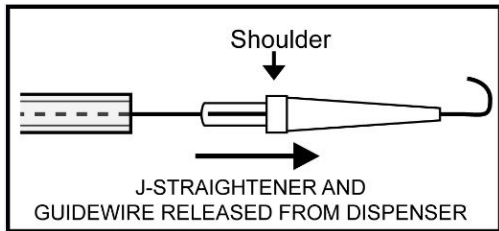

3. Holding the guidewire and the shank of the J-Straightener, fully withdraw the guidewire from the catheter.

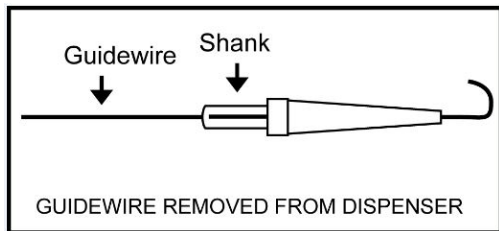

Note: Discard biohazard according to facility protocol.

CAUTION: DO NOT grasp and pull the guidewire prior to releasing the J-Straightener. Damage to the guidewire may occur if it is pulled against the restraint of the J-Straightener.

WARRANTY

MEDCOMP® WARRANTS THAT THIS PRODUCT WAS MANUFACTURED ACCORDING TO APPLICABLE STANDARDS AND SPECIFICATIONS. PATIENT CONDITION, CLINICAL TREATMENT, AND PRODUCT MAINTENANCE MAY EFFECT THE PERFORMANCE OF THIS PRODUCT. USE OF THIS PRODUCT SHOULD BE IN ACCORDANCE WITH THE INSTRUCTIONS PROVIDED AND AS DIRECTED BY THE PRESCRIBING PHYSICIAN.

Because of continuing product improvement, prices, specifications, and model availability are subject to change without notice. Medcomp® reserves the right to modify its products or contents without notice.

Medcomp® is a registered trademark of Medical Components, Inc.

Captive® J-Straightener is a registered trademark of Lake Region Manufacturing, Inc.

 **Medical Components, Inc.**

1499 Delp Drive

Harleysville, PA 19438 U.S.A.

Tel: 215-256-4201

Fax: 215-256-1787

www.medcompnet.com

EC REP

MPS Medical Product Service GmbH
Borngasse 20
35619 Braunfels
Germany

CE
2797