


DUO-JET ST

The unique *Duo-Jet ST* catheter overcomes several complications normally associated with standard "D" design double lumen acute access catheters. The fixed circular tube within a tube allows for a 360° arterial uptake and reduces the need for reversing blood lines. The design also increases catheter stability in-situ as compared to conventional double lumen catheters, which exhibit a "whipping" action during use.



DUO-JET ST ij

Much like the Duo-Jet ST, the *Duo-Jet ST ij* has a unique preformed curve design that overcomes patient discomfort, vessel trauma and catheter kinking when internal jugular access is indicated. This catheter is designed for a low, supraclavicular approach to the internal jugular vein. When inserted in this manner, the catheter will rest comfortably over the clavicle and the chest wall. Recent studies have demonstrated significantly reduced infection and bacteremia rates in internal jugular access with this catheter in the supraclavicular position.

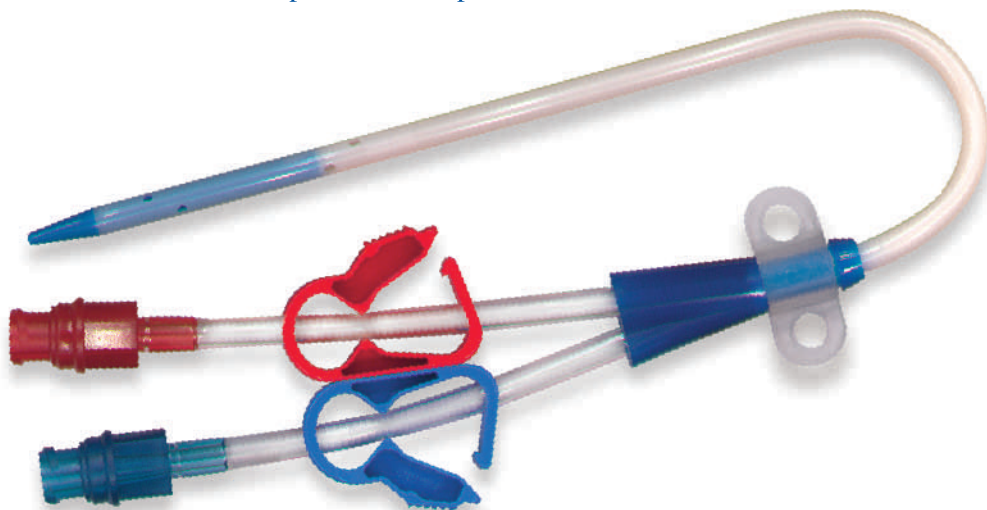

Features & Benefits

- Thermosensitive Flexi-Smooth® catheter material softens after insertion for increased patient safety.
- Soft, highly radiopaque tip reduces incidence of vessel damage and aids in confirming catheter tip replacement.
- Priming volume designation for ease of catheter maintenance.
- 360° venous ports decrease vessel turbulence while increasing blood flow.
- 360° arterial ports help eliminate vessel occlusion and arterial insufficiency.
- Seldinger Technique insertion may be applied.
- Rotatable suture wing provides secure anchoring.
- Durable high memory urethane extension lines will withstand repeated clamping.
- French size and length designations for ease of identification.
- Strain relief collar reduces kinking at the exit site.

Visit us on the Web at www.jetmedical.net

ORDERING INFORMATION

SETS 5/BOX	DESCRIPTION
Straight	
DJST114	12cm X 11.5FR Double Lumen Catheter
DJST116	15cm X 11.5FR Double Lumen Catheter
DJST118	20cm X 11.5FR Double Lumen Catheter
DJST119	24cm X 11.5FR Double Lumen Catheter
Curved (IJ)	
DJIJ114	12cm X 11.5FR Double Lumen IJ Catheter
DJIJ116	15cm X 11.5FR Double Lumen IJ Catheter
DJIJ118	20cm X 11.5FR Double Lumen IJ Catheter

Sets Contain:

- (1) Catheter
- (1) Introducer Needle
- (1) Vessel dilator
- (1) J/Flex Guidewire
- (2) Injection Caps

Jet Medical SA

Rue Fritz-Courvoisier 40
2300 La Chaux-de-Fonds
Switzerland
Tel: +41-32-9679021
Fax: +41-32-9679001
www.jetmedical.net
customerservice@jetmedical.net